



同田生物
TAUTO BIOTECH

地址：上海市张江高科技园区爱迪生路326号 邮编：201203

销售热线：400-820-3059

86-21-51320588-8008/8035/8026/8025


传真：86-21-51320502

QQ:800015916

MSN:Tauto_shanghai@hotmail.com

URL: www.tautobiotech.com

Mail: sales2@tautobiotech.com export2@tautobiotech.com

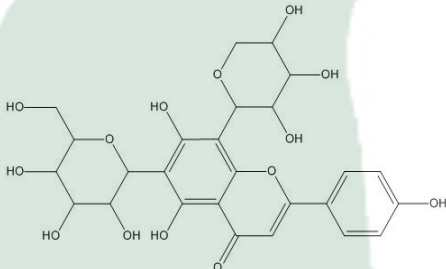PRODUCT ID	E-1182	Structure 
PRODUCT	Vicenin III/维采宁-3	
STANDARD TYPE	Quantitative Reference Standard/ 含量测定	
LOT NUMBER	15092122	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2015-09-21	
DATE OF REPORT	2015-09-21	

CHEMICAL NAME

6-C-Glucosyl-8-C-xylosylapigenin

OTHER NAME

CHEMICAL FORMULA

C₂₆H₂₈O₁₄

MOLECULAR WEIGHT(MW)

564.496

CAS NUMBER

59914-91-9

CHEMICAL FAMILY

FROM

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK298—2007	NLT ≥98.0%	Conforms
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK298—2007	NA	ND
Water	Q/TFFK298—2007	≤2%	Conforms
Appearance		Yellow Powder	Conforms

STORAGE CONDITIONS

STORAGE

Refer to the label

SHELF LIFE

2 Years


Analyst

Checker

Manager

Product Analysis Report

Product Name: Vicenin III **Chrom-Column:** C18
Lot Number: 15092122 **Mobile Phase:** MeOH: 0.1%H3PO4 (35: 65)
Detector: UV **Wavelength:** 270nm
Detect Time: 2015-09-21 **Inject Vol:** 20ul
Analyst: shuzhilin **Flow:** 1.0ml/min

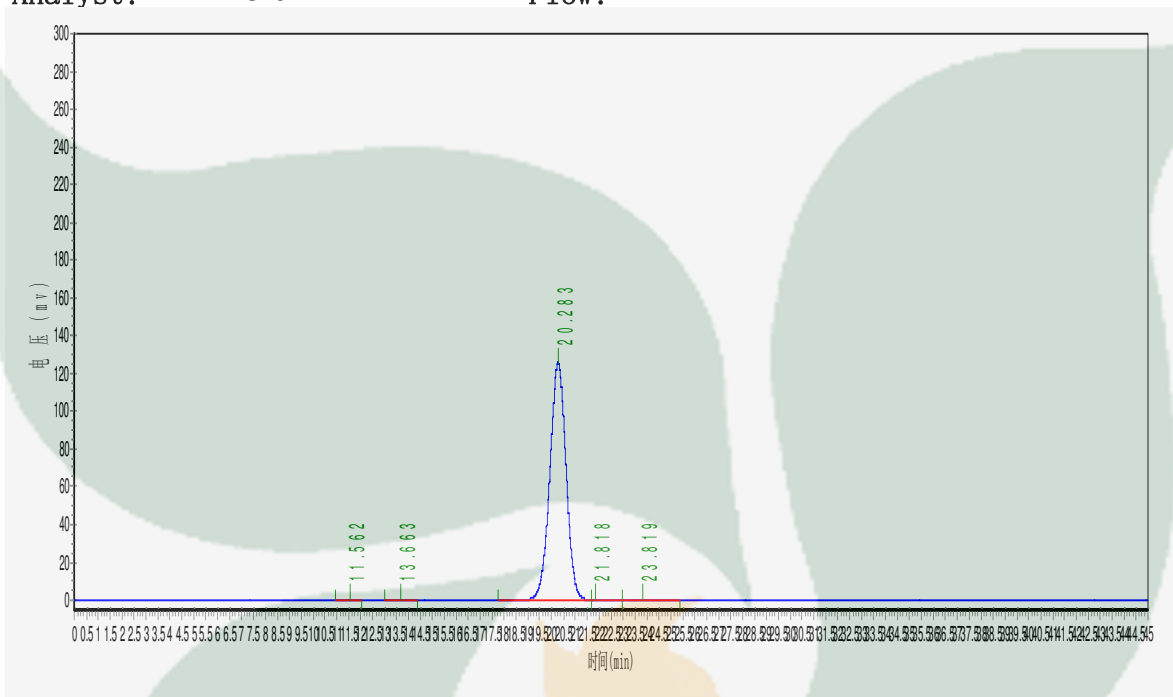

No.	Name	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1		11.562	451.877	13653.300	0.2271
2		13.663	268.046	9048.522	0.1505
3		20.283	125413.375	5953370.000	99.0035
4		21.818	422.537	17286.438	0.2875
5		23.819	348.837	19931.818	0.3315