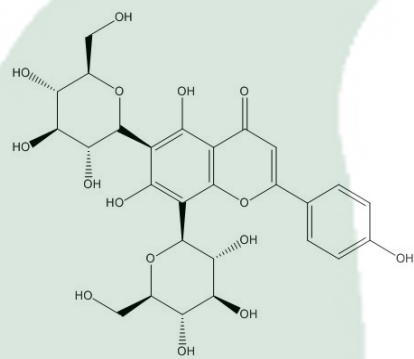


PRODUCT ID	E-1183	Structure 
PRODUCT	Vicenin-II/维采宁-2	
STANDARD TYPE	Quantitative Reference Material (含量测定)	
LOT NUMBER	16011122	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2016-01-11	
DATE OF REPORT	2016-01-11	

CHEMICAL NAME

OTHER NAME

CHEMICAL FORMULA

C₂₇H₃₀O₁₅

MOLECULAR WEIGHT(M)

594.519

CAS NUMBER

23666-13-9

CHEMICAL FAMILY

FROM

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK735—2010	NLT ≥98.0%	Conforms
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK735—2010	NA	ND
Water	Q/TFFK735—2010	≤2%	Conforms
Appearance		White Powder	Conforms

STORAGE CONDITIONS

STORAGE

Refer to the label

SHELF LIFE

2 Years


Analyst

Checker

Manager

Product Analysis Report

Product Name: Vicenin-II **Chrom-Column:** C18
Lot Number: 16011122 **Mobile Phase:** MeOH:0.1%H3PO4=35:65
Temperature: 25°C **Wavelength:** 272nm
Detect Time: 2016-01-11 **Inject Vol:** 20ul
Analyst: shuzhilin **Flow:** 1.0ml/min

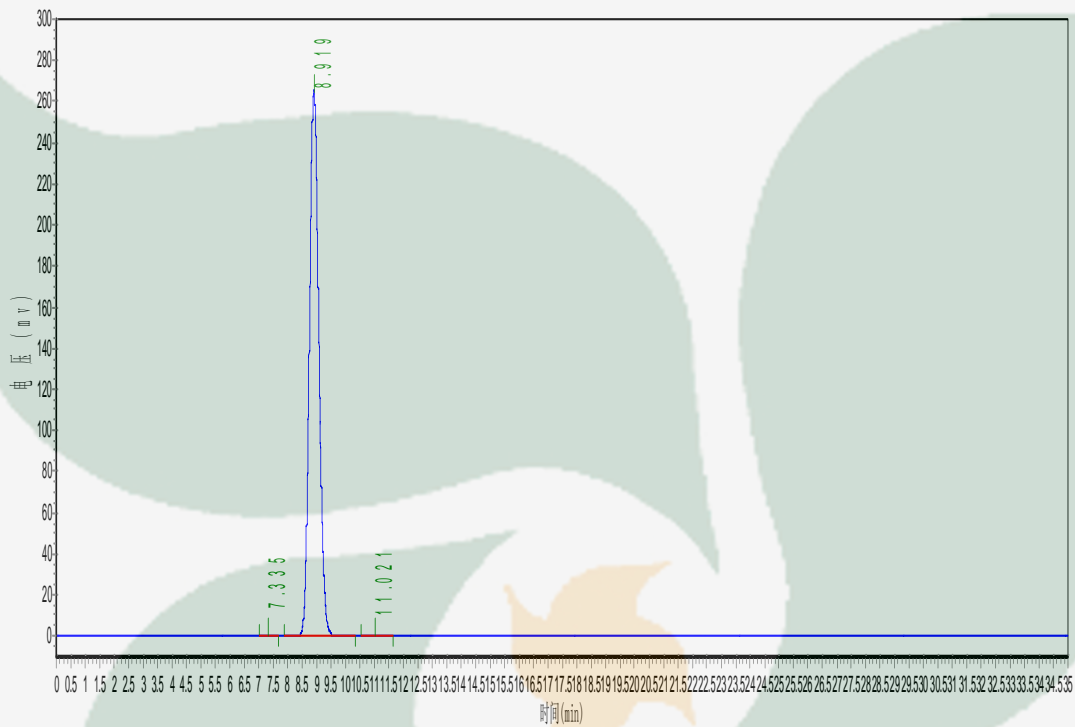

No.	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1	7.335	273.949	4946.100	0.0814
2	8.919	265199.656	6069280.500	99.8429
3	11.021	157.164	4600.800	0.0757