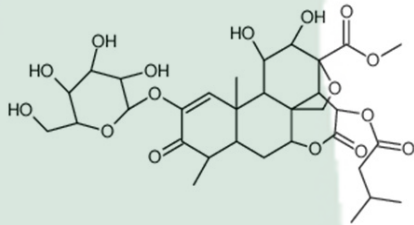


PRODUCT ID	E-2897	<p>Structure</p> 
PRODUCT	Yadanzioside A/鸦胆子内脂A	
STANDARD TYPE	Quantitative Reference Material (含量测定)	
LOT NUMBER	17062621	
ASSAY METHOD	HPLC,MS,NMR	
DATE OF SAMPLE	2017-06-26	
DATE OF REPORT	2017-06-26	

OTHER NAME

CHEMICAL FORMULA

C₃₂H₄₄O₁₆

MOLECULAR WEIGHT(MW)

684.263

CAS NUMBER

95258-15-4

CHEMICAL FAMILY

FROM

鸦胆子 *Brucea javanica*

ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK88—2007	NLT ≥98.0%	Conforms
NMR		Conforms	Conforms
MS		Conforms	Conforms
Residual Solvent	Q/TFFK88—2007	NA	ND
Water	Q/TFFK88—2007	≤2%	Conforms
Appearance		White Powder	Conforms

STORAGE CONDITIONS

STORAGE

Refer to the label

SHELF LIFE

2 Years


Analyst

Checker

Manager

Product Analysis Report

Product Name: Yadanzioside A **Chrom-Column:** C18
Lot Number: 17062621 **Mobile Phase:** ACN:H2O=10: 90
Detector: UV **Wavelength:** 210nm
Detect Time: 2017-06-26 **Inject Vol:** 20ul
Analyst: shuzhilin **Flow:** 1.0ml/min

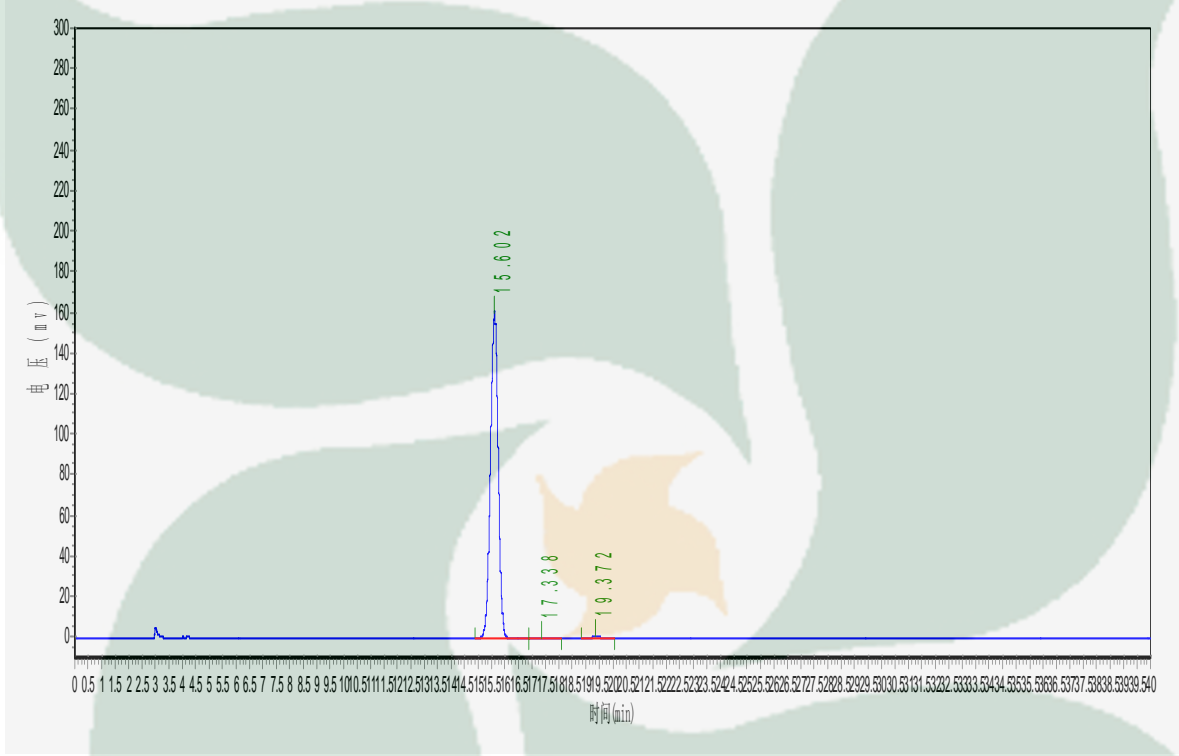

No.	Ret.Time min	Height mAU	Area mAU*min	Rel.Area %
1	15.602	160864.922	3242558.500	99.1759
2	17.338	143.884	4453.036	0.1362
3	19.372	966.866	22492.398	0.6879