



**同田生物**  
**TAUTO BIOTECH**

地址：上海市张江高科技园区爱迪生路326号 邮编：201203

销售热线：400-820-3059


86-21-51320588-8008/8035/8026/8025

传真：86-21-51320502

QQ:800015916

URL: www.tautobiotech.com

Mail: sales2@tautobiotech.com export2@tautobiotech.com
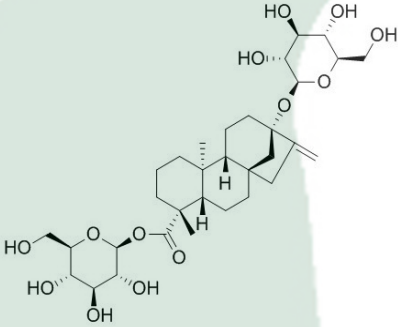
<b>PRODUCT ID</b>	E-1017	<b>Structure</b> 
<b>PRODUCT</b>	Rubusoside/甜叶悬钩子苷	
<b>STANDARD TYPE</b>	Quantitative Reference Standard/ 含量测定	
<b>LOT NUMBER</b>	18021431	
<b>ASSAY METHOD</b>	HPLC,MS,NMR	
<b>DATE OF SAMPLE</b>	2018-02-14	
<b>DATE OF REPORT</b>	2018-02-14	

<b>OTHER NAME</b>	NA
<b>CHEMICAL FORMULA</b>	C32H50O13
<b>MOLECULAR WEIGHT(MW)</b>	642.73
<b>CAS NUMBER</b>	64849-39-4
<b>CHEMICAL FAMILY</b>	Diterpenoids

#### ANALYTICAL RESULT

TEST	METHOD	SPECIFICATION	RESULT
HPLC	Q/TFFK95—2007	NLT ≥98.0%	98.6%
NMR		Conforms	Conforms
MS		Conforms	Conforms
<b>Residual Solvent</b>	Q/TFFK95—2007	NA	ND
<b>Water</b>	Q/TFFK95—2007	≤2%	Conforms
<b>Appearance</b>		White Powder	Conforms

#### STORAGE CONDITIONS

STORAGE Refer to the label  
SHELF LIFE 2 Years

Analyst

Checker

Manager

# Product Analysis Report

Product Name: Rubusoside

Chrom-Column: C18

Lot Number: 18021431

Mobile Phase: MeOH:0.1%H3PO4 (75:25)

Detector: UV


Wavelength: 208nm

Detect Time: 2018-02-14

Inject Vol: 20ul

Analyst: zhanghong

Flow: 1.0ml/min

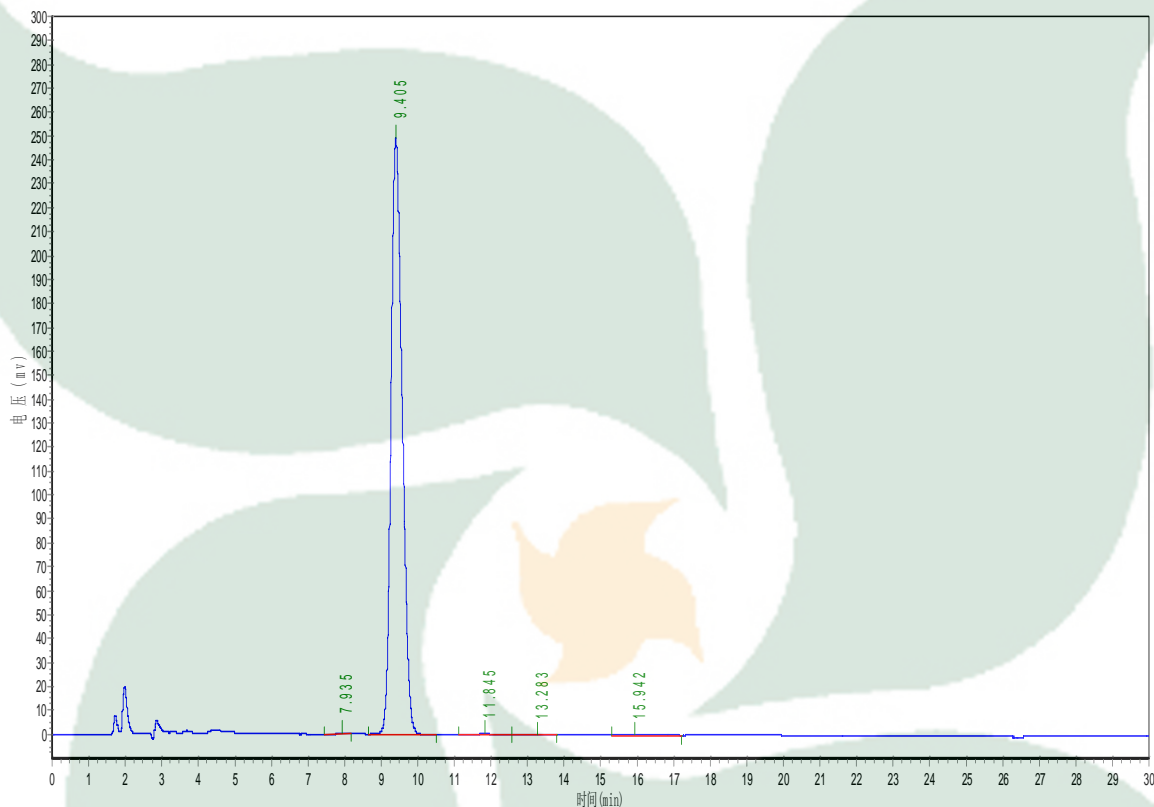

No.	Ret.Tim min	Height mAU	Area mAU*min	Rel.Area %
1	7.935	362.698	3674.650	0.0690
2	9.405	249519.516	5255697.000	98.6974
3	11.845	530.733	20702.900	0.3888
4	13.283	307.326	12249.147	0.2300
5	15.942	357.340	32734.729	0.6147